



Aesthetic, Plastic, Reconstructive Surgery and Dermatology

Specialized Operating Tables and Chairs

Welcome to BRUMABA

Reliable quality for your health

BRUMABA is a family-run company that specializes in the development and production of ergonomic operating tables, chairs and stools:

By developing high-quality, innovative products, we aim to improve the health and well-being of healthcare professionals and patients around the world.

The principle of our work is sustainability: BRUMABA products have been manufactured in our Bavarian homeland since 1980. The use of high-quality materials and reliable technology provide a proven long service life of the products.

Convince yourself of our noticeable advantage in ergonomics, comfort, reliability and efficiency.

Reliability that counts. Stability that lasts.
BRUMABA - your partner in quality.

MADE IN GERMANY



SINCE 1980



Operating chair
BALANCE ADVANCE

Patient recovery chair
CADDY



Operating table
VARIUS

Your advantages with BRUMABA

By Professionals, for Professionals

Their experience and procedures from everyday life in the operating theater are an important part of our product development. We have been working closely with surgeons and OR staff since 1980.

This is what makes BRUMABA so valuable for you.

- The products are specially designed to meet the requirements in Aesthetic, Plastic, Reconstructive Surgery and Dermatology.
- You and your team work with our operating tables and chairs in an ergonomic position at all times.
- Your patients are positioned professionally and comfortably during the entire surgical procedure - regardless of their size and weight.

Operating table
VARIUS



**Flexible to the
finest detail**



The art of surgery

The operating theater is the surgeon's stage

Every movement counts in the operating theater: Precision and concentration are crucial for the success of the operation and the well-being of the patient. Your surgical team works hand in hand, supported by essential instruments such as your BRUMABA operating table and operating chair. These help you to deliver your best performance at all times - and in the long term.

Thanks to the slim BRUMABA design, you enjoy optimum access to all areas of the patient while operating. This allows you to work closer to the operating field and adopt ergonomically ideal positions. Even areas, that are hard to reach, can be accessed effortlessly thanks to the wide range of adjustment options. This allows you to operate precisely in a comfortable position.

BRUMABA enables you to work ergonomically and perform precise operations. Experience maximum safety and efficiency for your everyday procedures with our solutions.



Ideal position for
facelift



**Maximum
flexibility:**
One operating
table for
all procedures



Hand surgery in
natural situation



Beach chair position
for breast surgery



Lateral support



Lithotomy position



Prone position



Operating table **VARIUS**

All in One - Our masterpiece

Standout features

- Two separately adjustable and spreadable leg section, for optimum access between the legs
- Airplane function: Tilt the patient towards you to work without straining your back
- Memory function: save up to 8 frequently used positions for efficient surgical procedures
- Premium, stable and durable stainless steel construction for highest quality "Made in Germany"
- Extra thick cushion for even better comfort during long procedures
- 1 - 2 weeks of battery operation per charge guarantees reliable OR operation
- Three-dimensionally adjustable head rest
- Extensive height adjustment with DT double telescopic column (625 - 1055 mm)
- Up to 300 kg patient weight with full functionality



Slim design
and legroom

Optional equipment

- Our recommendation: 4-EASY-DRIVE for easy adaptation to the room situation and cleaning
- Foot switch for sterile, intraoperative adjustment



Foot switch
for convenient
adjustment



Scan the QR code
and experience VARIUS live
in the product video






Airplane function
with 17°
Rotation angle



Ergonomic work
between the legs, e.g.
liposuction, intimate
surgery

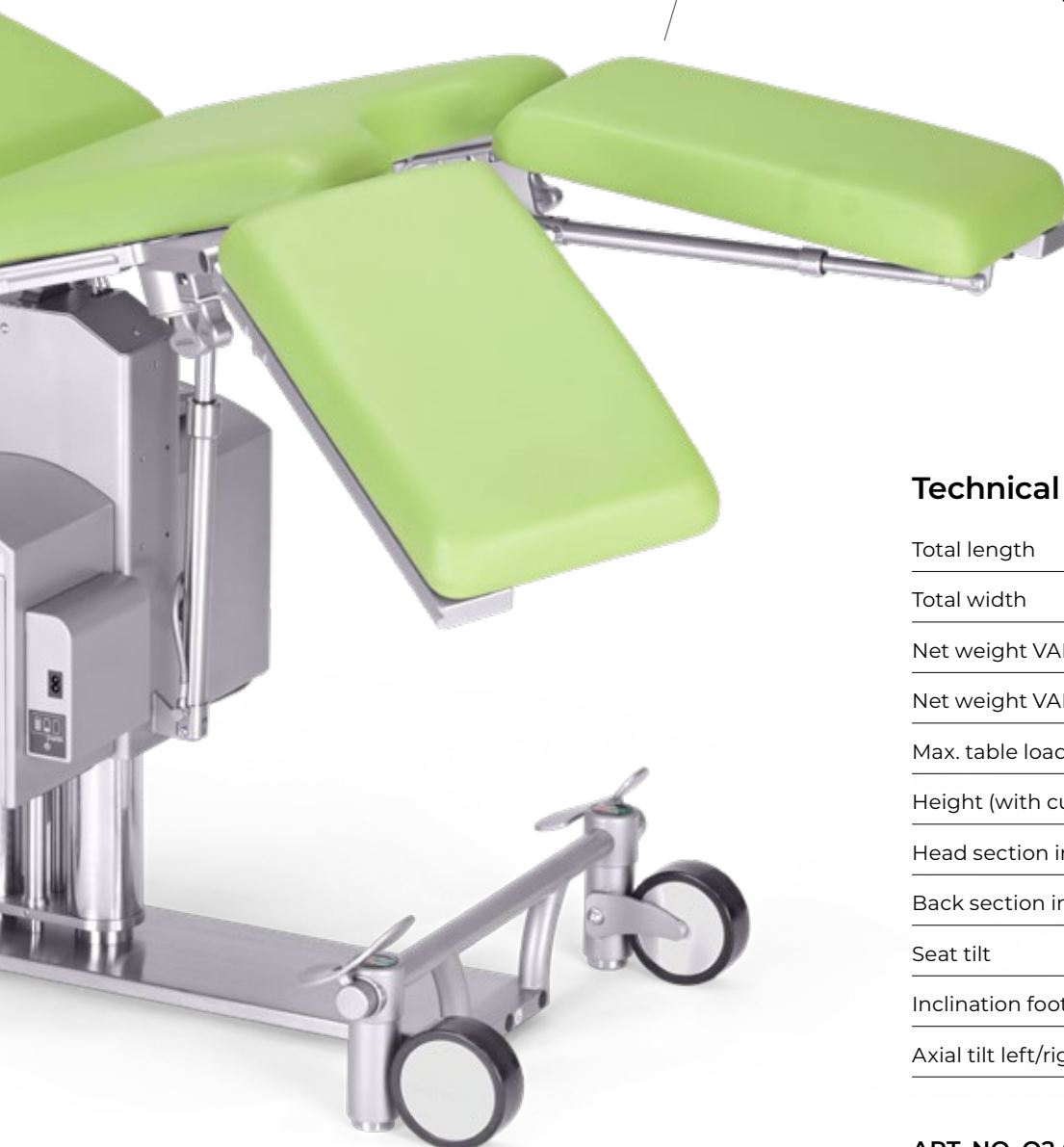

Technical data

Total length	1,770 - 1,940 mm
Total width	555 mm
Net weight VARIUS	205 kg
Net weight VARIUS N	195 kg
Max. table load	300 kg
Height (with cushion)	625 - 1055 mm
Head section inclination	31° to -44°
Back section inclination	61° to -14°
Seat tilt	43° to -31°
Inclination foot section	30° to -90°
Axial tilt left/right	17°

ART. NO. O2.02.31 VARIUS
ART. NO. O2.12.31 VARIUS N

Operating table **PRIMUS**

A unique operating experience

Standout features

- The airplane function enables perfect ergonomic positioning with longitudinal axis tilting
- Our three-dimensionally adjustable headrest elevates your surgical precision and comfort
- Crafted from high-quality stainless steel, ensuring stability and durability "Made in Germany"
- Flexibly adjustable armrests for relaxed positioning
- Extra thick cushions maximize patient comfort during long procedures
- Memory function allowing you to save up to 8 frequently used lying positions, improving operational efficiency
- Long-lasting battery life of 1-2 weeks, guaranteeing reliable performance
- Expanded height adjustment with DT double telescopic column (540 - 970 mm)
- Up to 300 kg / 660 lbs patient weight with full performance

Optional features

- Our recommendation: 2-GUIDE-DRIVE with two directional wheels for easy adjustment and cleaning in the operation room
- Foot switch feature for sterile and convenient intraoperative adjustments



Scan the QR code
and experience PRIMUS
live in the product video



Head shell for
e.g. rhinoplasty



Wireless battery
operation
for up to two weeks





High-quality stainless steel construction



Optimum ergonomics



Technical data

Total length	1,770 -1,940 mm
Overall width	555 mm
Deadweight	185 kg
Max. table load	300 kg
Height (incl. cushion)	540 - 970 mm
Inclination head section	31° to -44°
Inclination back section	72° to -17°
Inclination seat section	32° to -11°
Inclination foot section	30° to -90°

ART. NO. O2.01.31 PRIMUS
ART. NO. O2.11.31 PRIMUS N

OR table **GENIUS**

Excellent down to the smallest detail

Standout features

- Three-dimensionally adjustable head shell
- High-quality, stable and durable stainless steel construction for the highest quality "Made in Germany"
- Extra thick cushion for even better patient comfort during long procedures
- Memory function: storage of up to 8 frequently used positions for efficient surgical procedures
- 1 - 2 weeks of battery operation per charge guarantees reliable OR operation
- Height adjustment with ST telescopic column (590 - 880 mm)
- Up to 250 kg patient weight with full functionality

Optional equipment

- More extensive height adjustment with DT double telescopic column (530 - 960 mm)
- Foot switch for sterile, intraoperative adjustment
- Mounting rail on foot and head section for attaching accessories (ART. NO. Z1.151 Mounting rail 250 mm, Head/seat section)

ART. NO. G4.01.23 GENIUS DTA
ART. NO. G4.11.23 GENIUS DTN
ART. NO. G4.01.12 GENIUS STA
ART. NO. G4.11.12 GENIUS STN



Panel and rail construction




Wireless battery operation
for up to two weeks

Optional backrest
upholstery
with integrated armrests



The operating table shown contains special equipment.

Convenient storage

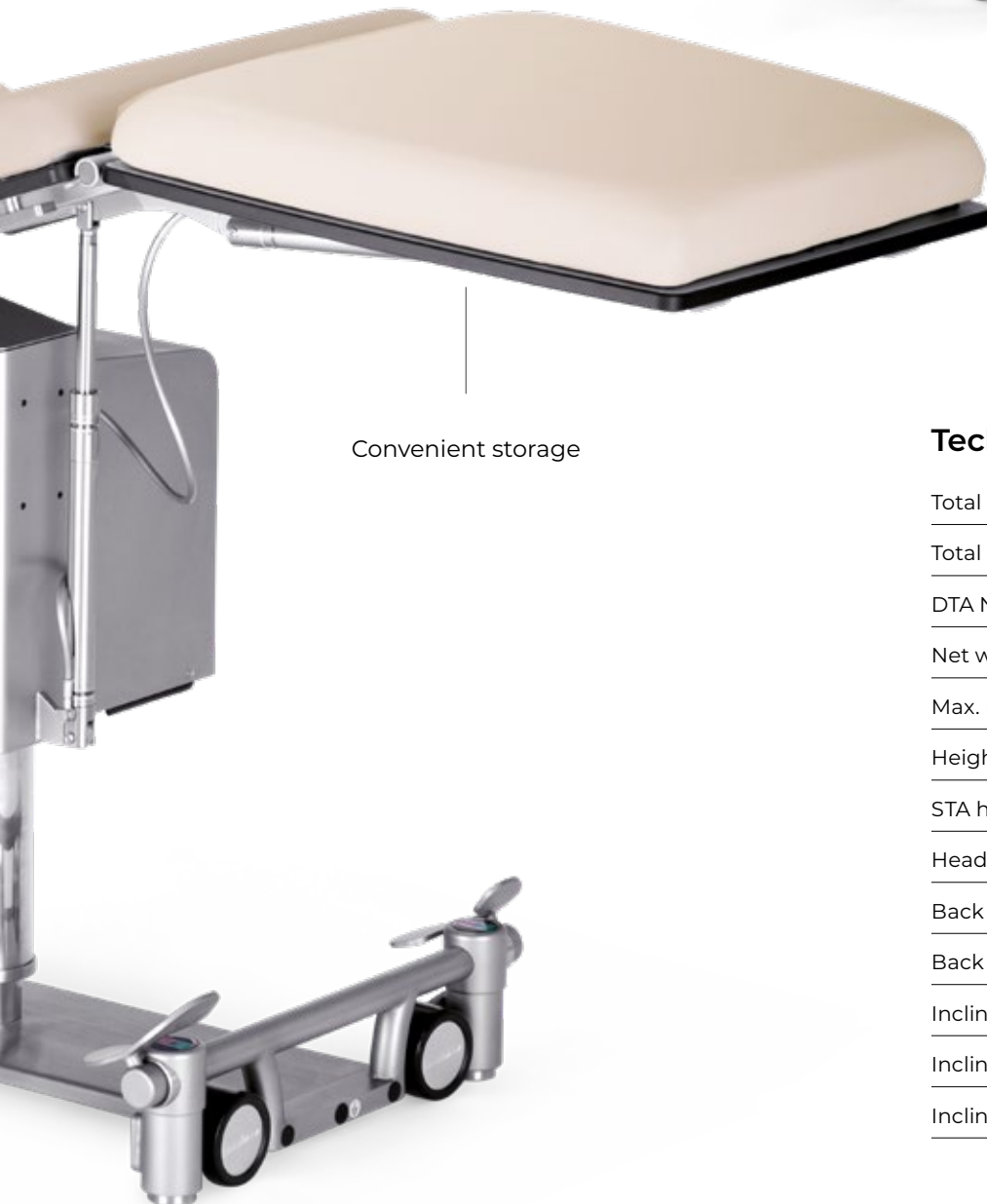

Technical data

Total length	1,803.5 - 1,991.5 mm
Total width	565 mm
DTA Net weight	160 kg
Net weight STA	140 kg
Max. Load capacity of the table	250 kg
Height DTA (incl. cushion)	530 - 960 mm
STA height (incl. cushion)	590 - 880 mm
Head section inclination	30° to -45°
Back section inclination ST	80° to -22°
Back section inclination DT	80° to -11°
Inclination seat section ST	39° to -9°
Inclination seat section DT	35 to -9°
Inclination foot section	28° to -90°

OR table **VENUS:**

2-in-1: treatment table and operating table

Standout features

- Magnetic neck pad and integrated nose slot for perfect tummy and back support
- High-quality, stable and durable stainless steel construction for the highest quality "Made in Germany"
- Extra thick cushion for even better patient comfort during long procedures
- Memory function: storage of up to 8 frequently used positions for efficient OR procedures
- 1 - 2 weeks of battery operation per charge guarantees reliable OR operation
- Height adjustment with ST telescopic column (630 - 910 mm)
- Up to 250 kg patient weight with full functionality



Optional equipment:

- More extensive height adjustment with DT double telescopic column (630 - 910 mm incl. cushion)
- Foot switch for sterile, intraoperative adjustment
- Mounting rail set for attaching accessories (ART. NO. Z1.304)

ART. NO. V1.01.12 VENUS STA

ART. NO. V1.01.23 VENUS DTA

ART. NO. V1.11.12 VENUS STN

ART. NO. V1.11.23 VENUS DTN



Prone position



Ergonomic positioning for surgery and treatment



Mounting rail set VENUS
(6 pcs.) for additional
accessories

ART.NO. Z1.304

Technical data

Total length	1,963 mm
Total width	660 mm
DTA Net weight	145 kg
Dead weight DTN	165 kg
Net weight STA	145 kg
Net weight STN	125 kg
Max. Load	250 kg
Height (with cushion)	630-1,060 mm
Height DTA (incl. cushion)	630-910 mm
Upholstery height	80 mm

Trolleys for operating tables

For maximum mobility and stability



Fixed base plate

Stable and Stationary: Ensuring the table remains securely in one place.

ART. NO. F1.101 PRIMUS

ART. NO. F1.103 GENIUS



2-GUIDE-DRIVE – GENIUS/PRIMUS

Easily move or lock the operating table with one lever. Features two directional and swivel castors, as well as a lift-brake system.

ART. NO. F1.203

OUR RECOMMENDATION



4-EASY-DRIVE

- 4 swivel wheels with connecting axle
- Lift-brake system for extremely smooth standing
- Intelligent directional guidance
- Maximum legroom
- Lightweight 120 mm wheels

ART. NO. F1.403



5-WHEEL DRIVE

- Selective 5th wheel for easy, precise movement in long corridors and around tight corners
- Lift-brake system for maximum stability
- Maximum legroom
- 120 mm easy-running wheels

ART. NO. F1.504

Controls and charging systems for operating tables



Foot switch type IV with retaining bracket

- Type IV foot switch with 8 memory positions

ART. NO. Z1.185



Foot control with holding bracket IV

For positioning the foot control on the chassis 4-EASY-DRIVE, 5-WHEEL-DRIVE.

ART. NO. Z1. 274



Manual control - VARIUS/PRIMUS

- Included as standard
- With 8 memory positions

ART. NO. Z1.118

**OUR
RECOMMENDATION**



Integrated charging system - PRIMUS/VARIUS

Charge your operating table once a week by plugging a cable directly into.

ART. NO. X.000269



Interchangeable battery system with external charging station - PRIMUS

Replacement battery for smooth non-stop operation, particularly suitable in large operating rooms with changing OR staff to ensure additional operational reliability.

ART. NO. Z1.134 Charging station

ART. NO. Z1.017 Battery 24 V
(one battery already included)



Rechargeable battery version 24 V – GENIUS

For independent operation of your operating table. Running time during normal operation of up to 2 weeks without recharging.

ART. NO. G4.0X.XX



Power supply 110-230 V

For stationary use, the mains operation is suitable for direct connection to a socket.

ART. NO. 02.11.31 PRIMUS N

ART. NO. G4.11.23 GENIUS DTN

ART. NO. G411.12 GENIUS STN

Accessories for operating tables

Head positioning



**Horseshoe-shaped
headrest, small**

200x210x80 mm

ART. NO. Z1.043A



Head plate, small

320x240 mm

ART. NO. Z1.109

Head plate, large

400x300 mm

ART. NO. Z1.003



**Hair transplant
support "neck support"**

ART. NO. Z1.032




Armrest OP PRO
with fixing strap and
length adjustment
ART. NO. Z1.331

Leg holder - OP,
with 2 ball-and-socket
joints (1 pair)
ART. NO. Z1.343



Hand surgery table
with telescopic support
ART. NO. Z1.351

Mounting rail 250mm,
solid – **GENIUS/VENUS**
Required for use of the
leg holder or the hand
surgery table with **GENIUS**
and **VENUS**
ART.NO. X.000087

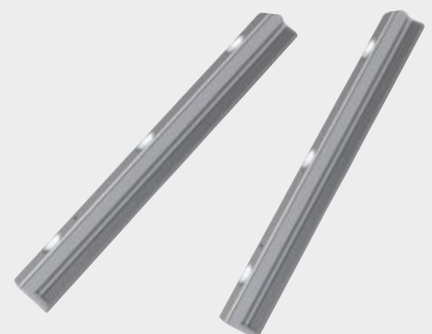

Accessories

For patient positioning like on cloud nine



**Back section -
extension, narrow**

Included with VARIUS, PRIMUS
and GENIUS



**Foot section extension /
side extension, narrow - OP**

ART. NO. Z1.267
PRIMUS/VARIUS

ART. NO. Z1.263
GENIUS



**Arm support lateral,
with fastening belt**

ART. NO. Z1.002



**Standing plate -
VARIUS/VENUS (1 pair)**

ART. NO. Z1.057

**Lateral positioning
support, deep**

180 x 300 x 50 mm

ART. NO. Z1.110



**Standing plate -
VARIUS/VENUS (1 pair)**

ART. NO. Z1.112




**Mattress viscoelastic,
short with
gynecological cutout**
ART. NO. Z1.213



**Mattress
viscoelastic**
ART. NO. Z1.171



Half roll cushion
ART. NO. Z1.100 (small)
ART. NO. Z1.101 (large)

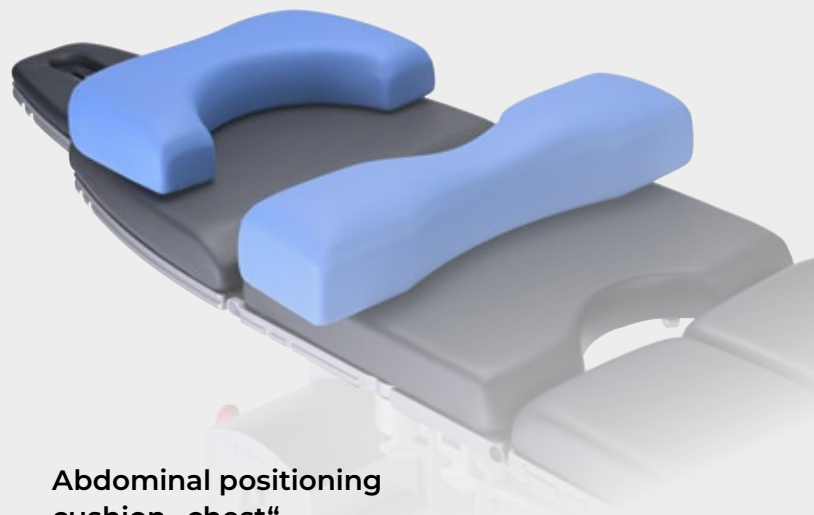

**Abdominal positioning
cushion „chest“**
ART. NO. Z1.345

**Abdominal positioning
cushion „hips“**
ART. NO. Z1.346



Lordosis support cushion
ART. NO. Z1.140



**Stomach positioning
cushion, trapezoidal
(1 piece)**
ART.NO. Z1.098

Instrument tables

Instrument table on the column

Plate size 600 x 400 mm


ART. NO. Z1.275

Plate size 700 x 500 mm

ART. NO. Z1.278



Infusion stand
on the column
ART. NO. Z1.280

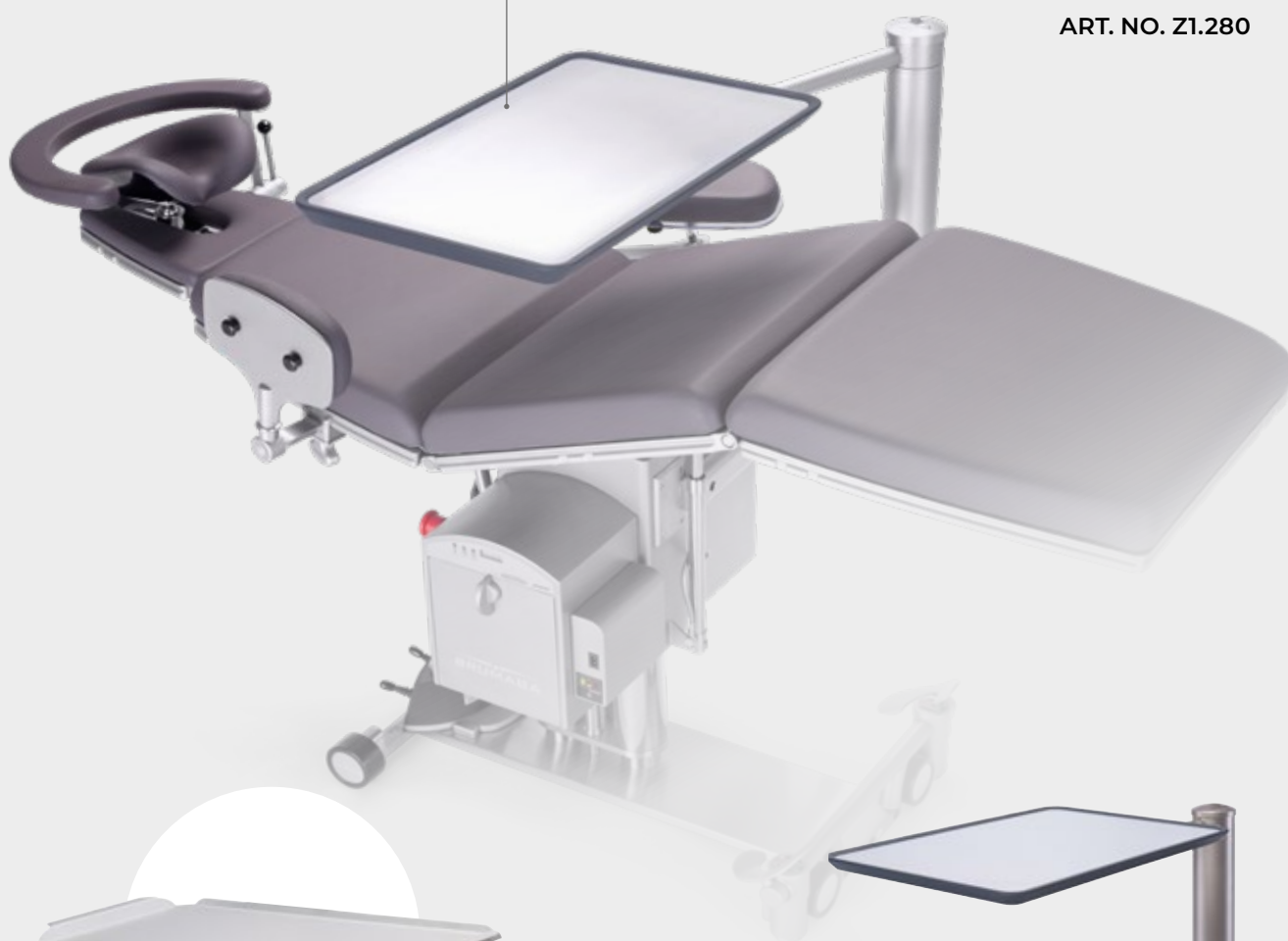

Aluminum shelf for
Instrument table
ART. NO. Z1.245

Mobile Instrument table

Plate size 600x400 mm

ART. NO. Z1.172

Plate size 700x500 mm

ART. NO. Z1.282



Remote control holders



Remote control holder
on
infusion stand
ART. NO. Z1.243



Remote control holder
FLEX, hinged
ART. NO. Z1.313




Remote control
holder, short
ART. NO. Z1.307



Magnetic plate for
remote control
ART. NO. Z1.045



Anesthesia accessories

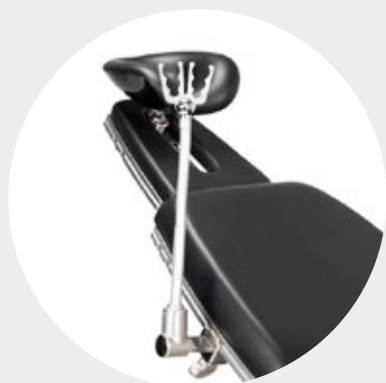

**Anesthetic tube holder
FLEX, foldable
ART. NO. Z1.334**



**Infusion stand on the
column, hinged
ART. NO. Z1.277**



**Infusion stand with
ball-and-socket joint
ART. NO. Z1.038**



**Anesthetic frame,
universal
ART. NO. Z1.015**

Protection, hygiene & assembly



Mounting bracket

ART. NO. Z1.004



**EasyCover
(150 pieces) disposable
table cover**

ART. NO. Z1.138



Paper roll holder

ART. NO. Z1.294 VARIUS

ART. NO. Z1.324 VENUS

ART. NO. Z1.051 PRIMUS/GENIUS



**Mounting
glider with
clamp**

ART. NO. Z1.093

**Protective Cover
for foot section**

ART. NO. Z1.332

ART. NO. Z1.148A
VARIUS



**Mounting adapter,
25x10 mm**

For attaching
accessories from other
manufacturers

ART. NO. Z1.047



**Mounting rail
– GENIUS/VENUS**

Included as standard
with GENIUS

250 mm Head and seat section
ART.NR. Z1.151

340 mm back and foot section
ART.NR. Z1.079

Treatment chair **ECURA**

Exclusive comfort meets stylish design

For over 45 years, BRUMABA has set the standard in surgical equipment. With ECURA we bring our full expertise to the treatment room – combining trusted BRUMABA quality with ergonomic design for fatigue-free working and exceptional client comfort.

Highlights

- Smart memory function: the handheld remote control can be used to save 4 frequently used reclining positions for maximum efficiency.
- Slim base for maximum legroom: the specially designed base allows greater freedom of movement and ensures a comfortable treatment position.
- The memory foam cushion, magnetic headrest, and ergonomic shape combine to provide a particularly comfortable and supportive reclining experience.
- Fast and quiet: The recliner can be quickly and quietly adjusted electrically in both height and position.

Technical Specifications

Min. height (incl. cushion)	638 mm
Max. height (incl. cushion)	889 mm
Table length	1,837 mm
Surface width	660 mm
Max. load capacity	200 kg

ART. NO. E1.00.00



Hand control



Convenient
access



Ideal
positioning



Ergonomic
working



Convenient
position



Foot control

Patient recovery chair **CADDY**

Design and comfort in a wake-up table

The CADDY is a 4-in-1 solution that functions as a chair, couch, transporter, and recovery chair. It is the perfect all-rounder that guarantees safe and comfortable patient positioning, coupled with easy transportation.

Standout features

- Robust frame crafted from high-quality stainless steel
- Two multifunctional handles on the upper edge of the backrest for ergonomic pushing and easy adjustment of the reclining surfaces
- Comfortable cushions
- Physiologically advantageous lying positions ensure stable circulation for your patients
- Easy transport thanks to particularly smooth-running wheels with steering brakes
- Fold-up armrests ensure convenient entry and better access to the patient
- Choice of light or dark wood

Optional accessories

- Armrest cushion with ball joint (ART. NR. Z1.070)
- Bodybelt, silicone (ART. NR. Z1.078)
- Caddy - Side protection (1 pair)
- Infusion stand (ART. NR. Z1.068)
- CADDY - wooden optic "beech light" standard (ART. NO. C1.XX1)
- CADDY - wooden optic "beech dark" (ART. NO. C1.XX1, surcharge)



Scan the QR code and experience CADDY live in the product video



Side protection

Armrests automatically adjusted to reclining position (foldable)

Lying surface infinitely adjustable



Technical data

Weight	67 kg
Max. dimensions (W x L x H) „Chair“	660 x 1100 x 1440 mm
Max. dimensions (W x L x H) „Lounge“	660 x 1900 x 600 mm
Seat dimensions (W x L)	460 x 490 mm
Seat height	560 mm
Max. load	180 kg

ART. NO. C1.02.13

Comfortable in any position



Chair position



Chaiselongue position



Shock position

BRUMABA OR stools and chairs

Invest in your health

BRUMABA operating stools and chairs can be flexibly adapted to your individual body dimensions and are the ideal complement to our operating tables.

Back protection minimizes physical strain during surgery and assisting at – a sustainable investment in your own health and that of your team.

The Standout features of BRUMABA operating chairs

- Optimum back support thanks to infinitely adjustable back section
- Precisely configurable armrests for perfect support during the operation
- Patented ergonomic seat surfaces for a consistently secure and comfortable fit
- Individually adjustable tilt angle of the seat for healthy bloodcirculation in your legs
- Adjustable seat height for an angle of $>90^\circ$ between your upper body and legs





Stool model
COMFORT



Backrest adjustment



Height adjustment with
foot release



Optimal
Back support

Ergonomic
Seat surface

Vein-friendly
Sitting

Operating chair
BALANCE

Microsurgery chair **BALANCE ADVANCE**

The comfortable and ergonomic solution for precision work

In close cooperation with microsurgeons, we have developed the BALANCE ADVANCE operating chair specifically for the requirements of microsurgery. The result: comfortable surgery with a perfect working position.





Experience unmatched comfort and ergonomics

- Three-directional adjustable armrests with lockable feature support shoulder and arm muscles
- Comfortable work position even during extended procedures
- Ergonomically shaped seat and back with patented seat inclination and high-quality cushions ensure a healthy sitting posture
- Easy height adjustment with foot release (gas spring)
- Comfortable chassis with antistatic wheels
- High-frequency seamless cushions and closed wheels for hygienic cleaning

BALANCE ADVANCE F
ART. NO. S3.10.13.5310

BALANCE ADVANCE F&B (with brake)
ART. NO. S3.11.13.5310

BALANCE ADVANCE SLIM F
ART. NO. S3.10.23.641.0

BALANCE ADVANCE SLIM F&B (with brake)
ART. NO. S3.12.23.641.0

Technical data

Min. height	555 mm
Max. height	745 mm
Chair cross diameter	610 mm
Seat cushion length	460 mm
Seat cushion width	420 mm

BALANCE ADVANCE SLIM

Min. height	490 mm
Max. height	630 mm
Seat cushion length	400 mm
Seat cushion width	360 mm

Operating chair **BALANCE**

Achieving perfect balance: Form, function, and ergonomics

BALANCE F

Standout features

- Body-adaptive seat and backrest for ergonomic support
- High-quality cushions for comfortable and healthy sitting
- The ergonomic seat design promotes healthy blood circulation in your legs, enabling you to sit comfortably for extended periods of time
- The cushions are seamlessly designed and the wheels are enclosed, allowing effortless cleaning and ensuring impeccable hygiene
- Optional column elevation (100 mm) (ART. NO. E.001558)



Height adjustment
with foot release,
ideal for sterile work

Technical data

Min. height	500 mm
Max. height	695 mm
Chair base diameter	610 mm
Seat cushion length	460 mm
Seat cushion width	420 mm

BALANCE F
ART. NO. S3.10.13.501.0

BALANCE F&B (with brake)
ART. NO. S3.11.13.501.0



Optional with brake

BALANCE SUPREME F

The BALANCE SUPREME is equipped with premium armrests that are padded and multifunctional.

- The armrests are adjustable in height, length, and can be moved to the side, providing a comfortable and customizable sitting experience
- Select the color of the armrest cushions to match the seat surface

BALANCE SUPREME F
ART. NO. S3.10.13.521.0

BALANCE SUPREME F & B (with brake)
ART. NO. S3.11.13.521.0

Technical data

Min. height	520 mm
Max. height	715 mm
Chair base diameter	610 mm
Seat cushion length	460 mm
Seat cushion width	420 mm

Operating stool **COMFORT**

The allrounder

Designed to enable medical staff to work with speed and agility, this maneuverable stool offers maximum flexibility for seamless use. With its practical foot-operated height adjustment, users can easily achieve the desired height while keeping their hands free.

Standout features

- Seamless seat cushions for enhanced comfort and easy cleaning, ensuring reliable hygiene
- Large, enclosed, antistatic wheels that facilitate fast and agile working, while maintaining optimal hygiene
- Optional locking brake for a secure and stable working position
- Optional 100 mm column elevation for increased versatility and customizable height adjustment (ART. NO. E.001558)



Technical data

Min. height	540 mm
Max. height	735 mm
Chair cross diameter	610 mm
Seat diameter	400 mm
Flexible height adjustment	510 - 710 mm



Height adjustment with foot release, ideal for sterile work



Available with COMFORT SWING seat cushion



COMFORT ASSIST H

With 360° swivelling backrest for ergonomic support in a variety of sitting positions.

Only available with height adjustment via manual release.

ART. NO. H3.10.21.103



Simple height adjustment of the backrest

Available COMFORT models

COMFORT F

Height adjustment by foot release

ART. NO. H3.10.11.100.0

COMFORT F/B

Height adjustment by foot release, with brake

ART. NO. H3.11.11.100.0



Optionally with brake

COMFORT SWING F

With COMFORT SWING seat cushion, height adjustment by foot release

ART. NO. H3.10.12.100.0

COMFORT SWING F/B

With COMFORT SWING seat cushion, height adjustment by foot release, with brake

ART. NO. H3.11.12.100.0



Height adjustment with hand release

COMFORT H

Height adjustment by hand release

ART. NO. H3.10.31.100.0

COMFORT SWING H

With COMFORT SWING seat cushion, height adjustment via hand release

ART. NO. H3.10.32.100.0

High-quality cushions

Maximize your comfort and hygiene

Viscoelastic comfort cushions

Experience the ultimate in pressure relief and comfort with our viscoelastic cushions, designed to effectively prevent bedsores. From start to finish, we prioritize your comfort and that of your patients, ensuring a seamless and relaxing procedure.




Elevate your hygiene standards to new heights

- The cushions can be easily removed for fast cleaning and disinfection
- Seamless and waterproof design to provide effective protection against bacteria
- High-quality durable artificial leather
- Our cushions are specially designed to meet all the hygiene requirements for operating rooms, ensuring the utmost safety and cleanliness for patients and staff
- Available for individual reordering, providing the ultimate flexibility for your seating needs

35+ Vibrant Cushion Colors

Unleash Your Creativity

Our cushion colors are available for operating tables, stools, chairs, and treatment tables.
This allows your BRUMABA products to integrate seamlessly into your facilities.

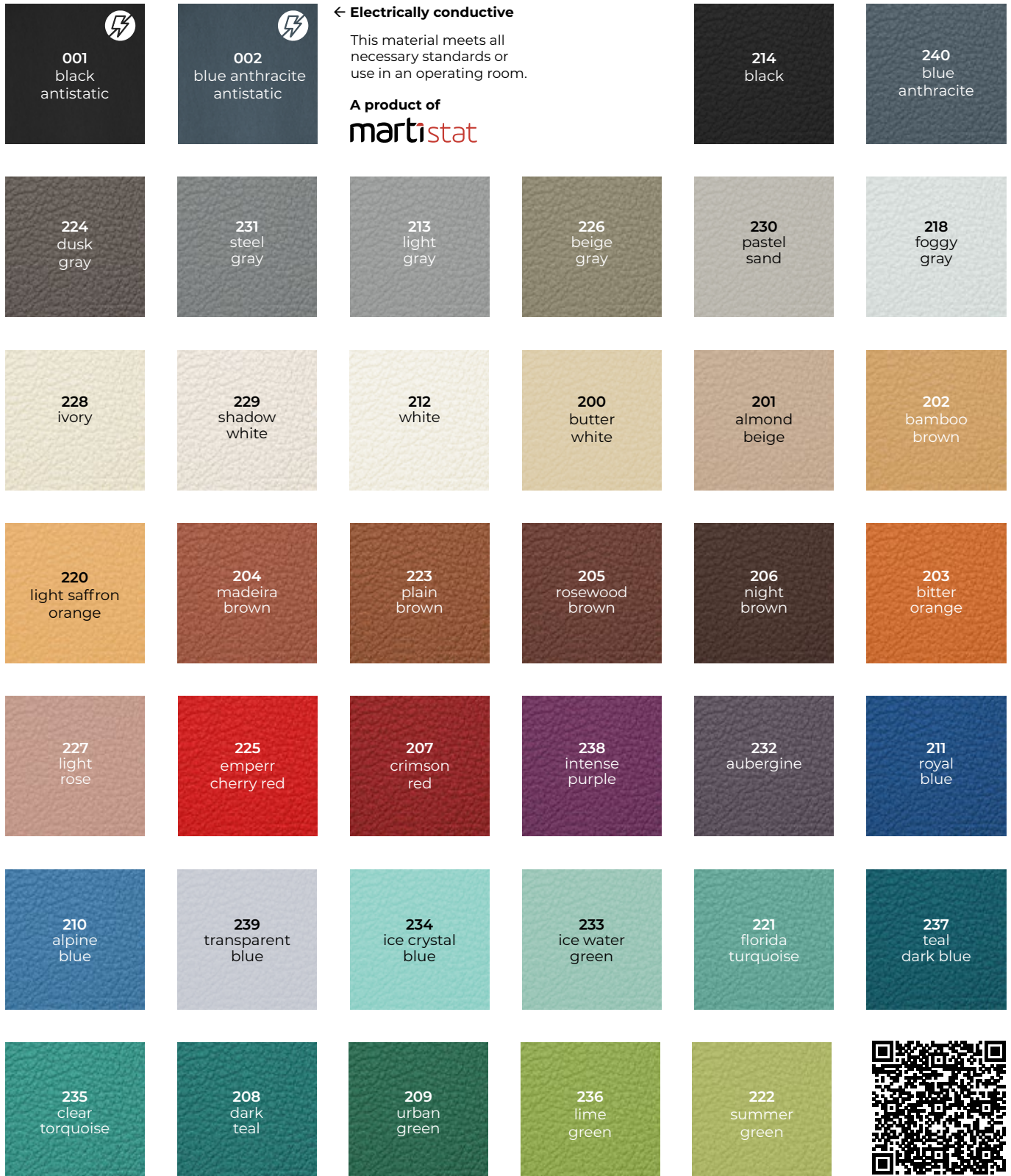

Scan QR code
for further information

Color swatches may differ from the original color tone.



BRUMABA GmbH

Bürgermeister-Graf-Ring 17
D-82538 Geretsried

Phone: +49 (0) 8171 / 2672 - 0

Fax: +49 (0) 8171 / 2672 - 10

Email: info@brumaba.de

Subject to change without notice.

© BRUMABA GmbH

KA.0015-2026-06-24-EN

www.brumaba.com

MADE IN GERMANY



SINCE 1980

BRUMABA
ERGONOMIC MEDICAL SOLUTIONS